TraceableLive®

Wi-Fi/Ethernet connected Traceable[®] Data Loggers with TraceableLIVE® Cloud Service

- Choose Wi-Fi or Ethernet connectivity
- New voice notification escalation*
- Apple Watch® integration
- New and improved user interface

*Requires a Premium subscription

Traceable is an Antylia Scientific company

0 1 1

191 3360-

4389° 1539, 1400,

75

Traceable

-

Storage Fridge

3.Fr MA PS: FO/ 1202, FO year, TA b

MA 45: 80/ 1202, 70 yest the be

noM anutsrag 100000018

lated At: May 07, 2021 (08:35 AM

....

3, Freeze

3, Freezer Roo

.....

0 160137592

Device 2CH 200450844

evices



Find Quality Products Online at:

www.GlobalTestSupply.com


TraceableLIVE[®] Cloud Service for remote monitoring and ease of mind

Refrigerator/Freezers • ULT Freezers • Chemical Storerooms • Greenhouses • Cryogenic (LN₂) Storage • Controlled Warehouses

Ideal for labs, research organizations, clinics, healthcare providers and anyone who requires that samples remain within set tolerances for temperature, humidity, barometric pressure, or CO₂. Protect your vaccines, pharmaceuticals, biological materials, or other product from degrading due to improper storage or transport conditions.

TraceableLIVE Cloud Service

- Cloud-based application—no local software required; web, Android[™], iOS data visibility
- Receive Apple Watch, Voice call*, mobile push, e-mail, text notifications
- Easy setup with mobile QR code
- Continuous data transmission and monitoring using SSL-encrypted SECURE connection
- Generate data reports in real time, data backup guaranteed
- View conditions and set alarm parameters on device remotely
- Expanded use with 21 CFR Part 11 compliance with premium subscription
- *Requires a Premium subscription

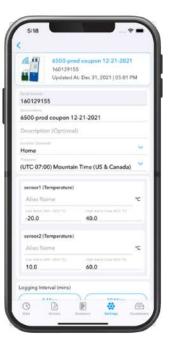



Reduce Your Risk

Real-time, remote monitoring—receive email, in app messages, or text alerts



Easy Mobile Setup Get connected to your devices and start monitoring in less than 2 minutes



Broad-Range Applications

Temperature, humidity, barometric pressure, CO₂

Find Quality Products Online at:


www.GlobalTestSupply.com

Personalize dashboards to immediately see data that is important to you.



evices (0, inclusive energy in Lander					
HIN DO	\$7 00		1285 L	G	bini 👘
Vaccine St. O juzzalian for	en tigan - 1		0	seation (the	atus. 1
time desc	- 3	Bator et			*
() General and a second second	. 0.	1000 (PAge 20.000 (20)	ne (skin és es)		9,
tob 01 @ Location: the	erronan. ; Oev	rice - CTP NCF	9	acation 10.12	cama 1
Laws Laws	-	4	Line		W hers Be not
to service a service service service	¥ .0=	anala-			ų,

Quickly find devices by name, location, or alarm status.

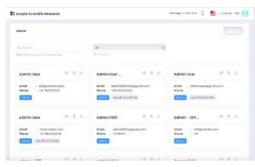

Easily manage user access and device details.

	Traceable live
Antyliz Scientific Research	antage surrene 🖉 🍯 Kapata dat
Ashish - 001 Potel - 01 - User Notification Settings	é sait see
ALLECATIONS	S Read Street Street Street
CANADA	Bins Bins Bins
London	المان هي المان هي المان هي المان الق
USA.	Directarian Carlos

You can even set notifications to escalate to other team members if alerts are not addressed.

Choose how you want to get notified of any alert status. Each location allows up to 3 notification types: email, in app messages, or text alerts.

2



At a glance see readings, alarm settings and statuses, Wi-Fi strength, battery level, and more.

View device data by channel, a specific time, or current time.

3

www.GlobalTestSupply.com

The sooner you know the better!

- Voice Notification Escalation now available for TraceableLIVE Premium subscribers
- TraceableLIVE Premium subscribers will now be able to set up an escalation path that generates an automated phone message for out-of-range conditions to mobile and landline phones
- Easy setup—just configure as an Escalation in the Escalations menu
- Ideal for high importance escalations



Get notified wherever you are it's all on your wrist!

- Accessing your TraceableLIVE devices couldn't be easier
- TraceableLIVE paid subscribers will be able to access their devices and notifications via Apple Watch



Find Quality Products Online at:

www.GlobalTestSupply.com

Meet compliance requirements with USB for VFC

- All TraceableLIVE users will be able to output secure PDF files as required by California code
- Download CSV files from all Traceable[®] Excursion-Trac[™] devices
- Upload USB data from Excursion-Trac data loggers



TraceableLive

Performance and speed improvements with app update

- App will automatically update on your phone OR you can manually update it if auto-updates are turned off
- All TraceableLIVE users will be able to receive the new updated app
- Update allows TraceableLIVE to better support users with future functionality improvements



5

Find Quality Products Online at:


www.GlobalTestSupply.com

TraceableLive

Prepaid Subscription Packages

**Meets 21 CFR Part 11 North American Compliance and Annex 11 for Europe.

Now you can get TraceableLIVE without monthly charges on your credit card. Simply order a subscription for 1 or 2 years. You can add it to your new TraceableLIVE product purchase or order it alone to add to one of your current TraceableLIVE devices.

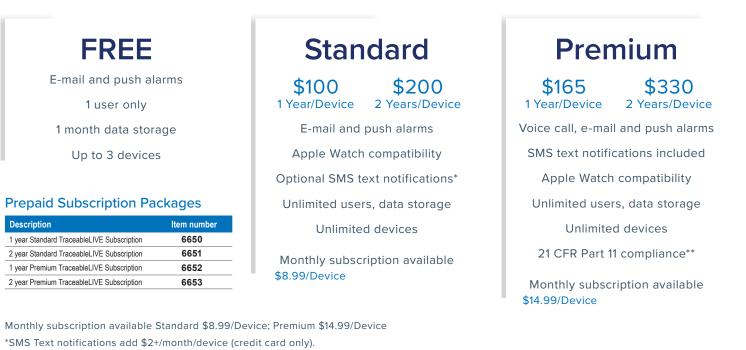

6

Find Quality Products Online at:

www.GlobalTestSupply.com

Traceable[®] Wi-Fi Data Loggers compatible with TraceableLIVE[®] Cloud Service^{*}

Get immediate alerts before you lose valuable samples

- Wi-Fi enabled digital data logger connects easily to network and provides continuous data transmission to the app interface
- Ideal for use in healthcare and vaccine monitoring (VFC compliant**)

What's Included: Probe with 10-ft (3-m) cable or USB dongle that is inserted into device, holder with magnetic back, four AAA alkaline batteries, and NIST-traceable calibration certificate with data.

Parameters	Range	Resolution	Accuracy	Probe	Item number
				2 glycol-filled temperature buffered bottles probes	6500
Temperature	–50.00 to 60.00°C (–58.00 to 140.00°F)	0.01°	±0.25°C	1 bottle & 1 bullet	6501
				2 plastic bullet waterproof probes	6502
				2 stainless steel 316 probes	6503
ULT	-90.00 to 105.00°C	0.01°	+0.1°C	1 platinum RTD sensor, stainless steel 316 probe	6510
Temperature	(-130.00 to 221.00°F)	0.01	±0.1 C	2 platinum RTD sensor, stainless steel 316 probes	6511
LN ₂	-200.00 to 105.00°C	0.01°	±0.25°C	1 platinum RTD sensor, stainless steel 316 probe	6512
Temperature	(-328.00 to 221.00°F)	0.01	±0.25 C	2 platinum RTD sensor, stainless steel 316 probes	6513
High	0.00 to 300.00°C	0.01°	±0.25°C	1 platinum RTD sensor, stainless steel 316 probe	6514
Temperature	(32.00 to 572.00°F)	0.01	±0.25 C	2 platinum RTD sensor, stainless steel 316 probes	6515
Temperature	0.00 to 55.00°C (32.00 to 131.00°F)	0.01°	±0.4°C	1 ambient dongle sensor	6520
Humidity	0.0 to 100%RH (noncondensing)	1%	±3% (20 to 80% RH), ±5% otherwise	1 external probe sensor on 3' cable	6521
CO ² Temperature Humidity	0 to 10,000 ppm 0.00 to 50.00°C (32.00 to 122.00°F) 0 to 95% (noncondensing)	1 ppm 0.01° 0.01% RH	4% reading + 100 ppm ± 0.4°C ±3% (5 to 75% RH), ±5% otherwise	1 ambient dongle sensor	6525
CO ² Temperature Humidity	0.00 to 20.00% 0.00 to 50.00°C (32.00 to 122.00°F) 0 to 95% (noncondensing)	0.01% 0.01° 0.01% RH	4% reading + 1% full scale ± 0.4°C ±3% (5 to 75% RH), ±5% otherwise	1 external probe sensor on 3' cable	6526
Barometric Pressure Temperature Humidity	500 to 1100 millibars (mb) 0 to 55°C (32 to 131°F) 0 to 100% RH (noncondensing)	1 mbar 0.01° 0.01% RH	±4 mbar ±0.4°C ±3% (5 to 75% RH), ±5% otherwise	1 ambient dongle sensor	6529



Traceable® Refrigerator/Freezer Data Logger Thermometer



Traceable® RH/Temperature Data Logger with Ambient sensor



Traceable® CO₂/RH/Temperature Data Logger with external probe 6526

*A TraceableLIVE account is required for TraceableLIVE services.

**Temperature-buffered bottle probes meet requirements for VFC (Vaccine for Children Program).

Find Quality Products Online at:

www.GlobalTestSupply.com

Traceable[®] Ethernet Data Logging Thermometers compatible with TraceableLIVE[®] Cloud Service^{*}

Always on, always connected

- Thermometer recalibration is accomplished by simply replacing the probe
- View current temperature conditions anywhere—cloud-based interface sends email, text, and push notifications
- Connect to your data via smartphone, tablet, or PC using TraceableLIVE cloud service*

These data logging thermometers with Ethernet offer an easy and secure way to connect your smartphone, tablet or personal computer remotely with your water baths, refrigerators, freezers and biological/pharmaceutical lab equipment.

Once setup is complete, you can securely monitor parameters 24/7, control your alarms remotely, view your data logging history, run reports in real time and give access to other team members. Each thermometer includes a replaceable, calibrated probe. Thermometer 99469-01 includes a bottle probe filled with nontoxic glycol protecting it from rapid temperature changes when the refrigerator or freezer door is opened. Thermometer 18004-37 comes with a stainless steel probe that measures a wide range of temperatures.

What's Included: Probe with 10-ft (3-m) cable, holder with magnetic back, Ethernet cable, universal power adapter and plug with 5-ft (1.5-m) cord, three AAA alkaline batteries, and NIST-traceable calibration certificate.

Description	Range	Accuracy	Item number
Ethernet Data Logging Thermometer with Bottle Probe		— ±0.25°C	6600
Replacement Bottle Probe (calibrated) (for 6600 only)	-50.00 to 60.00°C (-58.00 to 140.00°F)		6602
Ethernet Data Logging Thermometer with Stainless Steel Probe			6601
Replacement Stainless Steel Probe (calibrated) (for 6601 only)	90.00 to 105.00°C (-130.00 to 221.00°F)		6603

*A TraceableLIVE account is required for TraceableLIVE Cloud service.



Traceable® Ethernet Data Logging Refrigerator/Freezer Thermometer 6600



Top view



Find Quality Products Online at:

www.GlobalTestSupply.com